
User Manual

Model 9260

Dual Channel/Differential Power Amplifier Publication No. 041210

Rev 1.0

Tabor Electronics Ltd.

Tabor Electronics Ltd.
P.O. Box 404, Tel Hanan Israel 20302
Tel: +972-4-821-3393, FAX: +972-4-821-3388

PUBLICATION DATE: March 1, 2016

Copyright 2003 by Tabor Electronics Ltd. Printed in Israel. All rights reserved. This book or parts thereof may not be reproduced in any form without written permission of the publisher.

WARRANTY STATEMENT

Products sold by Tabor Electronics Ltd. are warranted to be free from defects in workmanship or materials. Tabor Electronics Ltd. will, at its option, either repair or replace any hardware products which prove to be defective during the warranty period. You are a valued customer. Our mission is to make any necessary repairs in a reliable and timely manner.

Duration of Warranty

The warranty period for this Tabor Electronics Ltd. hardware is three years, except software and firmware products designed for use with Tabor Electronics Ltd. Hardware is warranted not to fail to execute its programming instructions due to defect in materials or workmanship for a period of ninety (90) days from the date of delivery to the initial end user.

Return of Product

Authorization is required from Tabor Electronics before you send us your product for service or calibration. Call your nearest Tabor Electronics support facility. A list is located on the last page of this manual. If you are unsure where to call, contact Tabor Electronics Ltd. Tel Hanan, Israel at 972-4-821-3393 or via fax at 972-4-821-3388. We can be reached at: support@tabor.co.il

Limitation of Warranty

Tabor Electronics Ltd. shall be released from all obligations under this warranty in the event repairs or modifications are made by persons other than authorized Tabor Electronics service personnel or without the written consent of Tabor Electronics.

Tabor Electronics Ltd. expressly disclaims any liability to its customers, dealers and representatives and to users of its product, and to any other person or persons, for special or consequential damages of any kind and from any cause whatsoever arising out of or in any way connected with the manufacture, sale, handling, repair, maintenance, replacement or use of said products.

Representations and warranties made by any person including dealers and representatives of Tabor Electronics Ltd., which are inconsistent or in conflict with the terms of this warranty (including but not limited to the limitations of the liability of Tabor Electronics Ltd. as set forth above), shall not be binding upon Tabor Electronics Ltd. unless reduced to writing and approved by an officer of Tabor Electronics Ltd.

This document may contain flaws, omissions or typesetting errors. No warranty is granted nor liability assumed in relation thereto. The information contained herein is periodically updated and changes will be incorporated into subsequent editions. If you have encountered an error, please notify us at support@taborelec.com. All specifications are subject to change without prior notice.

Except as stated above, Tabor Electronics Ltd. makes no warranty, express or implied (either in fact or by operation of law), statutory or otherwise; and except to the extent stated above, Tabor Electronics Ltd. shall have no liability under any warranty, express or implied (either in fact or by operation of law), statutory or otherwise.

PROPRIETARY NOTICE

This document and the technical data herein disclosed, are proprietary to Tabor Electronics, and shall not, without express written permission of Tabor Electronics, be used, in whole or in part to solicit quotations from a competitive source or used for manufacture by anyone other than Tabor Electronics. The information herein has been developed at private expense, and may only be used for operation and maintenance reference purposes or for purposes of engineering evaluation and incorporation into technical specifications and other documents, which specify procurement of products from Tabor Electronics.

FOR YOUR SAFETY

Before undertaking any troubleshooting, maintenance or exploratory procedure, read carefully the **WARNINGS** and **CAUTION** notices.



CAUTION
RISK OF ELECTRICAL SHOCK
DO NOT OPEN



This equipment contains voltage hazardous to human life and safety, and is capable of inflicting personal injury.



If this instrument is to be powered from the AC line (mains) through an auto-transformer, ensure the common connector is connected to the neutral (earth pole) of the power supply.



Before operating the unit, ensure the conductor (green wire) is connected to the ground (earth) conductor of the power outlet. Do not use a two-conductor extension cord or a three-prong/two-prong adapter. This will defeat the protective feature of the third conductor in the power cord.



Maintenance and calibration procedures sometimes call for operation of the unit with power applied and protective covers removed. Read the procedures and heed warnings to avoid “live” circuit points.

Before operating this instrument:

1. Ensure the proper fuse is in place for the power source to operate.
2. Ensure all other devices connected to or in proximity to this instrument are properly grounded or connected to the protective third-wire earth ground.

If the instrument:

- fails to operate satisfactorily
- shows visible damage
- has been stored under unfavorable conditions
- has sustained stress

Do not operate until, performance is checked by qualified personnel.

DECLARATION OF CONFORMITY

We: Tabor Electronics Ltd.
9 Hatasia Street, Tel Hanan
ISRAEL 36888

declare, that the 45MHz 34Vp-p Dual / Differential Channel Signal Amplifier

Model 9260

meet the intent of Directive 89/336/EEC for Electromagnetic Compatibility and the requirements of the Low Voltage Directive 73/23/EEC amended by 93/68/EEC. Compliance was demonstrated to the following specifications as listed in the official Journal of the European Communities:

Safety:

IEC/EN 61010-1 2nd Edition:2001+ C1, C2

EMC:

EN 50081-1 Emissions:

EN 55022 - Radiated, Class B

EN 55022 - Conducted, Class B

EN 50082-1 Immunity:

IEC 801-2 (1991) - Electrostatic Discharge

IEC 801-3 / ENV50140 (1993) - RF Radiated

IEC 801-4 (1991) - Fast Transients

The tests were performed on a typical configuration.

Table of Contents

Chapter 1 - PORTRAYAL

What's in This Chapter	1-1
Introduction	1-1
9260 Feature Highlights	1-2
Conventions Used in this Manual	1-2
Functional	1-3
Description	1-3
Configuration Options	1-3
Specifications	1-4
Front Panel	1-4
Connectors	1-4
Front Panel	1-6
Indicators	1-6
Rear Panel Parts and Connectors	1-6
Operating	1-8
Instructions	1-8

Chapter 2 - INSTALLATION

Installation Overview	2-1
Unpacking and Initial Inspection	2-1
Safety Precautions	2-1
Performance Checks	2-3
Operating Environment	2-4
Power Requirements	2-4
Grounding Requirements	2-5
Long Term Storage or Repackaging for Shipment	2-5
Preparation for Use	2-6
Installation	2-6
Factory Configuration of the 9260 Settings	2-6
Field Configuration of the 9260 Settings	2-7

Chapter 3 - MAINTENANCE, PERFORMANCE CHECKS and ADJUSTMENTS

What's in This Chapter	3-1
Disassembly Instructions	3-1
Special Handling of Static Sensitive Devices	3-2
Cleaning	3-3
Repair and Replacement	3-3
Performance Checks.....	3-4
Environmental Conditions	3-4
Recommended Test Equipment	3-4
Gain Accuracy Tests	3-5
Amplifier Bandwidth	3-6
Amplitude Accuracy	3-6
Amplifier Pulse Response.....	3-7
Pulse Response Checks	3-7
Auxiliary Input.....	3-8
Auxiliary Offset Checks	3-8
Auxiliary Input Checks.....	3-9
Adjustment Procedures.....	3-10
Amplifier Gain Adjustments.....	3-10
Gain Adjustments.....	3-10
Pulse Response Adjustments	3-11
Pulse Adjustments	3-11

Appendix A - SPECIFICATIONS	A-1
--	------------

List of Figures

1-1, The Model 9260	1-3
1-2, The 9260 Rear Panel	1-6
2-1, Field Modification of Channel 1 Input Impedance	2-7

List of Tables

2-1, Input Impedance Jumpers.....	2-7
2-2, Input coupling Jumpers	2-8
2-3, Output coupling Jumpers	2-8
3-1, Recommended Test Equipment.....	3-5
3-2, Gain Accuracy Tests	3-6
3-3, Output Bandwidth Tests	3-7
3-4, Square Wave Characteristics.....	3-8
3-5, Square Wave Characteristics.....	3-9

This page intentionally left blank

Chapter 1

PORTRAYAL

What's in This Chapter

This chapter contains general and functional description of the Model 9260 dual/differential wideband power amplifier. It also describes the front and rear panel connectors and its operational modes. It also and provides description of all features available with the instrument.

Introduction

The Model 9260, as shown in Figure 1-1, is a bench-top, 2U, half 19" rack size, fully metal case dual channel amplifier. The instrument can be configured to be used as two, single-ended and independent channels, or as a one input with two differential outputs.

The inputs to the amplifiers be configured to match different source impedances such as 50 Ω , 75 Ω , or 1 M Ω . The output is driven from a 0 Ω source.

There are three inputs for each channel:

1. Main input. This input is located on the front panel and is normally used for signal inputs.
2. Auxiliary input. This input is located on the rear panel and can be used as a summing input.
3. DC Offset input. This input is also located on the rear panel and can be used for offsetting the signal level within the specified output level window.

The outputs are located on the front panel. There are two outputs, one for each channel. When the 9260 is configured as two separate amplifiers, the outputs generate amplified signals within the range of 34 Vp-p.



Note

The output of the 9260 tends to oscillate if left open. The same may occur if the signal is applied to a high impedance source.

The bandwidth of the amplifier is around 30 MHz for large signals. Small signal bandwidth can reach 45 MHz.

The 9260 can be configured as a differential amplifier. In this case, the channel 2 input is disabled and channel 1 input is amplified and

distributed differentially to both outputs. In this case, channel 1 output generates in-phase signal while channel 2 outputs an inverted signal that has exactly 180° phase offset to the normal output. Full amplitude and bandwidth is preserved when the 9260 operates in differential mode. Using the differential mode, the 9260 does not sacrifice accuracy, nor does it sacrifice bandwidth.

The 9260 has two additional inputs for each channel allowing summation of two signals and providing an external control of DC level offset. These inputs are accessible from the rear panel only.

9260 Feature Highlights

- Bench-top, all metal case
- Large signal bandwidth to 30 MHz
- Small signal bandwidth to 45 MHz
- High amplitude to 34Vp-p (high impedance)
- Peak current to 1A
- Continuous output power to 10 W
- Low distortion

Conventions Used in this Manual

The following symbols may appear in this manual:



Note

A Note contains information relating to the use of this product



TIP

A Tip contains information relating to the performance of this product



Caution

A Caution contains information that should be followed to avoid personal damage to the instrument or the equipment connected to it.



Warning

A Warning alerts you to a potential hazard. Failure to adhere to the statement in a WARNING message could result in personal injury.



Figure 1-1, The Model 9260

Functional Description

Detailed functional description of the features, operation and options available with the 9260 is given in the following paragraphs. The wideband amplifier can be ordered with different configurations such as input impedance, gain, etc. Therefore read the following description carefully and make sure your amplifier is configured correctly for your application before you connect the cables to your circuits.

Configuration Options

The Model 9260 must be ordered from the factory already configured for your application. This manual has no schematics and no instructions how to modify the amplifier for other configurations as any configuration change, without full engineering supervision, may affect the performance of the amplifier. Below, you'll find a list of optional configurations for the amplifier:

Basic Configuration – defines if the 9260 will operate as two single-ended and independent amplifiers, or as single input with differential outputs.

Input Impedance – determines the matching source impedance at the input connectors. Three options are available: 50Ω, 75Ω and 1MΩ. Input impedance can be configured for the front panel main inputs and for the rear panel auxiliary inputs.

Gain – specifies gain magnitude of the input signal. Factory default setting is 10. Custom gain can also be specified. Note that some characteristics of the output section may change for gain setting above 10.

As explained above, all options must be specified at the time of your purchase and the 9260 is supplied fully configured. Reconfiguration of fielded instruments can be done by qualified and trained persons only.

Specifications

Instrument specifications are listed in Appendix A. These specifications are the performance standards or limits against which the instrument is tested. Specifications apply under the following conditions: output terminated into load impedance, after 30 minutes of warm up time, and within a temperature range of 20°C to 30°C. Specifications outside of the temperature range are degraded by 0.1% per °C.

Front Panel Connectors

The 9260 has 4 BNC connectors on its front panel, two for each channel. Two are marked INPUT and the other two are marked OUTPUT. These connectors are described below.

Channel 1 Input

The input connector accepts signals within the range of DC to over 35 MHz and amplifies them by a fixed gain. The gain is normally x10 however, this number may be different if you ordered the amplifier with different gain setting. Input impedance is factory pre-conditioned and has one of the following values: 50Ω, or 1MΩ. 1MΩ input impedance can be used for low frequency signals (up to 100kHz) however, higher frequencies require 50Ω termination at the input of the amplifier to eliminate standing waves and reflections in the input cable, which can cause excessive ringing and aberrations at the output.

The amplifier input cannot tolerate high voltage on its 50Ω resistance. Therefore, before applying the cable to the input connector, make sure your signal will not exceed input rating, as specified in Appendix A of this manual.



Note

Channel 1 Input only is used when the Model 9260 is configured to operate in differential mode.

Channel 2 Input

This input connector accepts signals within the range of DC to over 35 MHz and amplifies them by a fixed gain. The gain is normally x10 however, this number may be different if you ordered the amplifier with different gain setting. Input impedance is factory pre-conditioned and has one of the following values: 50 Ω , or 1M Ω . 1M Ω input impedance can be used for low frequency signals (up to 100kHz) however, higher frequencies require 50 Ω termination at the input of the amplifier to eliminate standing waves and reflections in the input cable, which can cause excessive ringing and aberrations at the output.

The amplifier input cannot tolerate high voltage on its 50 Ω resistance. Therefore, before applying the cable to the input connector, make sure your signal will not exceed input rating, as specified in Appendix A of this manual.



Note

Channel 2 Input is disabled when the 9260 is configured to operate in differential output mode.

Channel 1 Output

The channel 1 output connector outputs amplified signals. When the 9260 is configured as two separate amplifiers, this output connector generates amplified signals that are applied to the channel 1 input connector. If the instrument is configured for differential output, this channel generates normal, non-inverted signals while the other output connector generates the inverted signal at 180° phase offset.

Gain at this output is fixed at 10 (or another gain factor that was specified at the time of your purchase). Output source impedance is 0 Ω .



WARNING

Applying the output signal on inductive or capacitive loads may damage the amplifier.

Channel 2 Output

The channel 2 output connector outputs amplified signals. When the 9260 is configured as two separate amplifiers, this output generates amplified signals that are applied to the channel 2 input connector. If the instrument is configured for differential output, this channel generates an inverted signal at 180° phase offset while channel 1 output generates the normal, non-inverted.

Gain at this output is fixed at 10 (or another gain factor that was specified at the time of your purchase). Output source impedance is 0 Ω .

**WARNING**

Applying the output signal on inductive or capacitive loads may damage the amplifier.

Front Panel Indicators

The 9260 has 2 indicators on its front panel. The POWER LED will light as soon as you press the switch to power up the 9260. An illuminated power light designates power is applied to the instrument and you should expect to have signal at its output connector(s).

An LED at the center of the front panel is marked Differential ON. When this light is on, channel 2 input is disabled and signal applied to the channel 1 input is routed to both channel 1 and channel 2 outputs, except the signal is now differential. The differential mode is selected with a rear panel switch.

Rear Panel Parts and Connectors

The rear panel is rarely used for normal operation however, there are some connections you must do before you can operate the amplifier for example, connecting the mains power to the instrument. There are other connectors and parts on the rear panel; These are discussed in the following paragraphs. Refer to Figure 1-2 throughout the following description.



Figure 1-2, The 9260 Rear Panel

Line Receptacle and Fuse

Power is connected to the 9260 through the line receptacle. The amplifier accepts any voltage from 80 to 265Vac and there is no need to select the voltage range between different countries. Instructions how to connect the line cord and how to replace the line fuse is given in Chapter 2. If a fuse blows, make sure you replace it with the same type and rating to avoid possible damage to the product from unsuitable fuse value.

Auxiliary Inputs Channel 1

There are two Channel 1 auxiliary inputs on the rear panel marked INPUT and OFFSET.

The INPUT connector can be used as summing input to the front panel signal. It also can be used as the main input if connectors are available on the back only. However, the output is always from the front panel connectors. Observe the input limitations as specified in Appendix A to avoid damage to the amplifier circuit.

The OFFSET connector is used for dc signals only. Use this input to offset your signal to any direction, positive or negative however, always observe the positive and negative rail limitations as specified in Appendix A to avoid damaging the amplifier circuits.

With normal configuration, this input is used for driving the front panel output for channel 1 only. When the 9260 is configured for differential outputs, channel 2 auxiliary inputs are disabled and only channel 1 signal is amplified.

Auxiliary Inputs Channel 2

There are two Channel 2 auxiliary inputs on the rear panel marked INPUT and OFFSET.

The INPUT connector can be used as summing input to the front panel signal. It also can be used as the main input if connectors are available on the back only. However, the output is always from the front panel connectors. Observe the input limitations as specified in Appendix A to avoid damage to the amplifier circuit.

The OFFSET connector is used for dc signals only. Use this input to offset your signal to any direction, positive or negative however, always observe the positive and negative rail limitations as specified in Appendix A to avoid damaging the amplifier circuits.

With normal configuration, this input is used for driving the front panel output for channel 2 only. When the 9260 is configured for differential outputs, these inputs are disabled and only channel 1 signal is amplified.

Differential Outputs Switch

The Differential Outputs switch has two actions. In normal mode, it separates between the two channels so each channel can be used separately. In Differential Mode, this switch disables channel 2 inputs and converts the internal circuits to output differential signals. Channel 1 is used for outputting the normal output and channel 2 is used for outputting the inverted signal. When this switch is depressed, the front panel led illuminates, designating that the differential mode has been selected.

Grounding Considerations

Understanding how to connect your ground path could be critical to preserving the integrity of your output signal. If you are using a single-ended output then it will probably be safe for you to connect the circuit ground to case ground. Always bear in mind the following warning:



WARNING

Input and output grounds are tied together to case ground and therefore, it is absolutely forbidden to connect the output ground to a different level than the input ground. Failure to adhere to this limitation may damage the 9260 and the surrounding equipment connected to its I/O connectors.

Operating Instructions

Being a passive device, there are no controls, nor computer programming required to operate the Model 9260. The following procedure is recommended for proper operation of the power amplifier:

1. Make sure your box is configured for input and output impedance and gain
2. Follow the installation instructions given in Chapter 2 of this manual
3. Connect the output terminal to your load
4. Connect the input terminal to your source
5. Turn on power to your 9260



WARNING

There is no switch control to turn 9260 amplification on and off and therefore, the amplifier is active immediately after you power it up. Always make sure your load is protected from inadvertent power up conditions before you turn on your 9260.

Chapter 2

CONFIGURING the INSTRUMENT

Installation Overview

This chapter contains information and instructions necessary to prepare the Model 9260 for operation. Details are provided for initial inspection, grounding safety requirements, repackaging instructions for storage or shipment, installation information and Ethernet address configuration.

Unpacking and Initial Inspection

Unpacking and handling of the generator requires normal precautions and procedures applicable to handling of sensitive electronic equipment. The contents of all shipping containers should be checked for included accessories and certified against the packing slip to determine that the shipment is complete.

Safety Precautions

The following safety precautions should be observed before using this product. Although some instruments and accessories would normally be used with non-hazardous voltages, there are situations where hazardous conditions may be present.



CAUTION

This product is intended for use by qualified persons who recognize shock hazards and are familiar with the safety precautions required to avoid possible injury. Read the operating information carefully before using the product.

Exercise extreme caution when a shock hazard is present. Lethal voltage may be present on power cables, connector jacks, or test fixtures. The American National Standard Institute (ANSI) states that a shock hazard exists when voltage levels greater than 30V RMS, 42.4V peak or 60 VDC are present.



WARNING

For maximum safety, do not touch the product, test cables, or any other instrument parts while power is applied to the circuit under test. **ALWAYS** remove power from the entire test system before connecting cables or jumpers, installing or removing cards from the computer, or making internal changes such as changing the module address.

Do not touch any object that could provide a current path to the common side of the circuit under test or power line (earth) ground. Always keep your hands dry while handling the instrument.

When using test fixtures, keep the lid closed while power is applied to the device under test. Carefully read the Safety Precautions instructions that are supplied with your test fixtures.

Before performing any maintenance, disconnect the line cord and all test cables. Only qualified service personnel should perform maintenance.

Performance Checks

The instrument has been inspected for mechanical and electrical performance before shipment from the factory. It is free of physical defects and in perfect electrical order. Check the instrument for damage in transit and perform the electrical procedures outlined in the section entitled **Unpacking and Initial Inspection**.

Operating Environment

The 9260 is intended for operation on the bench or inside a rack of instruments. It is intended for indoor use only and should be operated in a clean, dry environment with an ambient temperature within the range of 0 °C to 50 °C.



WARNING

The 9260 must not be operated in explosive, dusty, or wet atmospheres. Avoid installation of the module close to strong magnetic fields.

The design of the 9260 has been verified to conform to EN 61010-1 2nd addition safety standard per the following limits: Installation (Overvoltage) Category I (Measuring terminals) Pollution Degree 2.

Installation (Overvoltage) Category I refers to signal level, which is applicable for equipment measuring terminals that are connected to source circuits in which measures are taken to limit transient voltages to an appropriately low level.

Pollution Degree 2 refers to an operating environment where normally only dry non-conductive pollution occurs. Occasionally a temporary conductivity caused by condensation must be expected.

Power Requirements

The amplifier may be operated from a wide range of mains voltage 90 to 264Vac. Voltage selection is automatic and does not require switch setting. The instrument operates over the power mains frequency range of 48 to 63Hz. Always verify that the operating power mains voltage is the same as that specified on the rear panel.

The 9260 should be operated from a power source with its neutral at or near ground (earth potential). The instrument is not intended for operation from two phases of a multi-phase ac system or across the legs of a single-phase, three-wire ac power system. Crest factor (ratio of peak voltage to rms.) should be typically within the range of 1.3 to 1.6 at 10% of the nominal rms. mains voltage.

Grounding Requirements

To ensure the safety of operating personnel, the U.S. O.S.H.A. (Occupational Safety and Health) requirement and good engineering practice mandate that the instrument panel and enclosure be “earth” grounded. Although BNC housings are isolated from the front panel, the metal part is connected to earth ground.

**WARNING**

Do not attempt to float the output from ground as it may damage the Model 9260 and your equipment.

**WARNING**

Input and output grounds are tied together and therefore, it is absolutely forbidden to connect the output ground to a different level than the input ground. Failure to adhere to this limitation may damage the 9260 and the surrounding equipment connected to its I/O connectors.

Long Term Storage or Repackaging for Shipment

If the instrument is to be stored for a long period of time or shipped to a service center, proceed as directed below. If repacking procedures are not clear to you or, if you have questions, contact your nearest Tabor Electronics Representative, or the Tabor Electronics Customer Service Department.

1. Repack the instrument using the wrappings, packing material and accessories originally shipped with the unit. If the original container is not available, purchase replacement materials.
2. Be sure the carton is well sealed with strong tape or metal straps.
3. Mark the carton with the model and serial number. If it is to be shipped, show sending and return address on two sides of the box.

**NOTE**

If the instrument is to be shipped to Tabor Electronics for calibration or repair, attach a tag to the instrument identifying the owner. Note the problem, symptoms, and service or repair desired. Record the model and serial number of the instrument. Show the RMA (Returned Materials Authorization) order as well as the date and method of shipment. ALWAYS OBTAIN AN RMA NUMBER FROM THE FACTORY BEFORE SHIPPING THE 9260 TO TABOR ELECTRONICS.

Preparation for Use Preparation for use includes removing the instrument from the container box and connecting the cables to its input and output connectors.

Installation

If this instrument is intended to be installed in a rack, it must be installed in a way that clears air passage to its cooling fans. For inspection and normal bench operation, place the instrument on the bench in such a way that will clear any obstructions to its rear fan to ensure proper airflow.



CAUTION

Using the 9260 without proper airflow will result in damage to the instrument.

Factory Configuration of the 9260 Settings

When you order the Model 9260, you should provide details how you want this product configured. There is some amount of flexibility in the configuration before the instrument is shipped from the factory. The following are factory configured settings:

Gain

The default gain setting is x10.

Input Impedance

The default input impedance is 50 Ω . The 9260 can be ordered with different input impedance values such as 75 Ω and 1M Ω . The 1M Ω is recommended for low frequency operation only as the high input impedance degrades the frequency response of the amplifier. Appendix A specifies the bandwidth for the various source and input impedances.

Source (output) Impedance

The default source (output) impedance is 0 Ω .

Input/Output Coupling

The default input/output coupling is DC. The 9260 can be ordered AC coupled. DC coupling is recommended for signals that have dc components with low frequencies operation. Appendix A specifies the frequency limitation for inputs/outputs configured to AC coupling.

Single-ended/Differential Outputs

The default configuration is single-ended. With this configuration, there are two independent channels that amplify signals that are applied to the channel 1 and 2 inputs. The 9260 can be ordered with differential outputs; In this case, channel 2 input is disabled and channel 1 output is amplified and routed differentially to both channels 1 and 2. Channel 1 outputs the normal signal and channel 2 outputs the inverted phase signal. Appendix A specifies the properties of the 9260 when configured for differential outputs operation.

Field Configuration of the 9260 Settings

Limited flexibility is available for you, as the user, to change 9260 configuration in the field. Changing configuration involves the opening of the top cover and therefore, it is suggested that this operation be completely avoided or, if absolutely necessary, should be carried out by factory-trained person.

Use the following instructions to re-configure 9260 settings in the field:

Remove the Top Cover

There are four screws bolting the top cover to the case. These are located on both sides of the cover. Identify and remove these screws. Use Philips screwdriver only to remove the screws. Hold the cover from both sides and pull upwards until the top is clear of the case sides. Once the top cover is open, you'll be able to access jumpers that set the different configuration.

Configure the Input Impedance

Before you change the 9260 input impedance settings, identify first the jumper location for all of the inputs. Figure 2-1 shows an example of the channel 1 input impedance jumpers. Place the jumpers as required according to the following table:

Table 2-1, Input Impedance Jumpers

Impedance	CH1	CH2	Aux CH1	Aux CH2
50Ω	J23	J5	J29	J14
75Ω	J24	J11	J28	J13
1MΩ	J22	J6	J30	J15

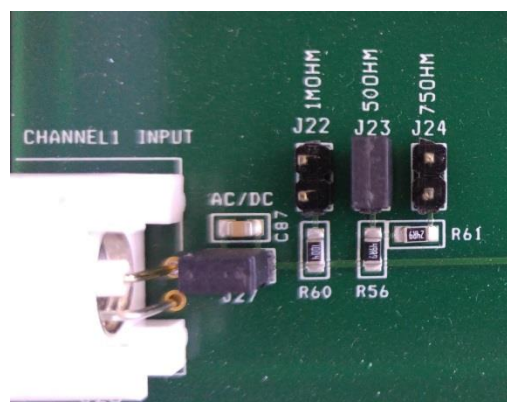


Figure 2-1, Field Modification of Channel 1 Input Impedance

Configure Input/Output Coupling

The default input/output coupling is DC. Before you change the 9260 input/output coupling settings, identify first the jumper location for all of the inputs and outputs. When the jumpers are on the links, the path is DC coupled. Remove the jumpers from the link to modify the settings to AC. Place or remove the jumpers as required according to table 2-3 and 2-4.

Table 2-2, Input coupling Jumpers

Input	CH1	CH2	Aux CH1	Aux CH2
Link	J27	J4	J31	J16

Table 2-3, Output coupling Jumpers

Impedance	CH1	CH2
Link	J19	J8

Configure Single-ended or Differential Outputs

The default configuration is single-ended. If you order the amplifier with differential outputs, channel 2 input is disabled and channel 1 output is amplified and routed differentially to both channels 1 and 2. Channel 1 outputs the normal signal and channel 2 outputs the inverted phase signal.

The rear panel, as shown in Figure 1-2 has a push-push switch. Pressing the switch once, will activate the differential mode and will cause a light on the rear panel to illuminate. Pressing the switch again will disable the differential mode and allow each amplifier to operate separately.

This page was intentionally left blank

Chapter 3

MAINTENANCE, PERFORMANCE CHECKS and ADJUSTMENTS

What's in This Chapter

This chapter provides maintenance and service information, performance tests, and the procedures necessary to adjust and troubleshoot the Model 9260 Universal Waveform Generator.



WARNING

The procedures described in this section are for use only by qualified service personnel. Many of the steps covered in this section may expose the individual to potentially lethal voltages that could result in personal injury or death if normal safety precautions are not observed.



CAUTION

ALWAYS PERFORM DISASSEMBLY, REPAIR AND CLEANING AT A STATIC SAFE WORKSTATION.

Disassembly Instructions

If it is necessary to troubleshoot the instrument or replace a component, use the following procedure to remove the side panels:

1. Using a Phillips head screwdriver, remove the screws from the top and bottom covers.
2. Carefully lift the top cover from its back end and slide the cover to the rear to clear the front panel spring latch. Do the same for the bottom. After removing the covers from the instrument, access the component side for calibration and checks, and the solder side when replacing components.
3. When replacing the top and bottom covers, reverse the above procedure.

Special Handling of Static Sensitive Devices

CMOS devices are designed to operate at very high impedance levels for low power consumption. As a result, any normal static charge that builds up on your person or clothing may be sufficient to destroy these devices if they are not handled properly. When handling such devices, use the precautions described below to avoid damaging them:

1. CMOS IC's should be transported and handled only in containers specially designed to prevent static build-up. Typically, these parts are received in static-protected containers of plastic or foam. Keep these devices in their original containers until ready for installation.
2. Ground yourself with a suitable wrist strap. Remove the devices from the protective containers only at a properly grounded workstation.
3. Remove a device by grasping the body; do not touch the pins.
4. Any printed circuit board into which the device is to be inserted must also be grounded to the bench or table.
5. Use only anti-static type solder sucker.
6. Use only grounded soldering irons.

Once the device is installed on the PC board, the device is adequately protected and normal handling may resume.

Cleaning

The Model 9260 should be cleaned as often as operating conditions require. To clean the instrument, use the following procedure:

1. Thoroughly clean the inside and outside of the instrument.
2. When cleaning inaccessible areas, remove dust with low-pressure compressed air or a vacuum cleaner.
3. Use alcohol applied with a cleaning brush to remove accumulation of dirt or grease from connector contacts and component terminals.
4. Clean the exterior of the instrument and the front panel with a mild detergent mixed with water, applying the solution with a soft, lint-free cloth.

Repair and Replacement

Repair and replacement of electrical and mechanical parts must be accomplished with great care. Printed circuit boards can become warped, cracked or burnt from excessive heat or mechanical stress. The following repair techniques are suggested to avoid inadvertent destruction or degradation of parts and assemblies:

1. Use a 60/40 solder and temperature-controlled 35 - 40 watt pencil-type soldering iron on the circuit board. The tip of the iron should be clean and properly tinned for best heat transfer to the solder joint. A higher wattage soldering iron may separate the circuit from the base material.
2. Keep the soldering iron in contact with the PC board for a minimum time to avoid damage to the components or printed conductors.
3. To de-solder components, use a commercial "solder sucker" or a solder-removing SOLDER - WICK, size 3.
4. Always replace a component with an exact duplicate as specified in the parts list.

Performance Checks

The performance of the 9260 should be checked to verify proper operation of the instrument and should normally be used:

1. As a part of the incoming inspection of the instrument specifications;
2. As part of the troubleshooting procedure;
3. After any repair or adjustment before returning the instrument to regular service.

Environmental Conditions

Tests should be performed under laboratory conditions having an ambient temperature of 25°C, $\pm 5^\circ\text{C}$ and at relative humidity of less than 80%. If the instrument has been subjected to conditions outside these ranges, allow at least one additional hour for the instrument to stabilize before beginning the adjustment procedure.

Warm-up Period

Most instruments are subject to small amount of drifts when first turned on. To ensure accuracy, turn on the power to the Model 9260 and allow it to warm-up for at least 10 minutes before beginning the performance test procedure.

9260 Configuration for the Tests

The 9260 is tested using the factory recommended configuration. The instrument should be configured as follows for the performance tests: Input impedance = 50 Ω ; Output impedance = 50 Ω ; Gain = 10; Coupling = DC. Other configurations can be tested but slight degradation of performance should be considered, as specified in Appendix A

Recommended Test Equipment

Recommended test equipment for troubleshooting, calibration and performance checking is listed below. Test instruments other than those listed may be used only if their specifications equal or exceed the required characteristics.

Table 3-1, Recommended Test Equipment

Equipment	Model No.	Manufacturer
Oscilloscope	DSO9404A	Keysight
Digital Multimeter	2000	Keithley
Spectrum Analyzer	FSP3	R&S
Pulse Generator	PM8571A	Tabor Electronics

Performance Check Procedures

Use the following procedures to check the Model 9260 against the specifications. A complete set of specifications is listed in Appendix A. The following paragraphs show how to set up the instrument for the test, what the specifications for the tested function are, and what acceptable limits for the test are. If the instrument fails to perform within the specified limits, the instrument must be calibrated or tested to find the source of the problem.

Gain Accuracy Tests

Equipment: DMM, PM8571A Pulse Waveform Generator

Preparation:

1. Configure the DMM as follows:

Function: ACV

Range: As required by the test

2. Connect the 9260 input to the PM8571A output
3. Connect 9260 Channel output to the DMM input via 50Ω/2W Termination
4. Configure the PM8571A as follows:
 - Frequency: 10 kHz
 - Outputs: On
 - Amplitude: As required for the test

Test Procedure

1. Perform Gain Accuracy tests using Table 3-2

Table 3-2, Gain Accuracy Tests

PM8571A Amplitude Setting	Error Limits	DMM Reading	Pass	Fail
0.5 Vp-p	1.767 V, ±71 mV			
1 Vp-p	3.535 V, ±141 mV			
2 Vp-p	7.071 V, ±282 mV			

Amplifier Bandwidth

Amplifier bandwidth checks the bandwidth of the output.

Amplitude Accuracy

Equipment: Oscilloscope, PM8571A Pulse Waveform Generator

Preparation:

1. Configure the Oscilloscope as follows:
 - Time Base: 200 μ s/div
 - Amplitude: 5 V/div
2. Connect the 9260 input to the PM8571A output
3. Connect 9260 Channel outputs to the Oscilloscope input via 50 Ω /2W Termination
4. Configure the PM8571A as follows:
 - Amplitude: 2 V
 - Outputs: On
 - Frequency: As specified in Table 3-3

Test Procedure

1. Using the variable vertical adjustment on the oscilloscope, adjust the vertical trace to show exactly 6 vertical divisions
2. Perform bandwidth tests using Table 3-3

Table 3-3, Output Bandwidth Tests

PM8571A Frequency Setting	Error Limits	Oscilloscope Reading	Pass	Fail
1 MHz	6 Divisions			
10 MHz	6 \pm 0.5 Divisions			
20 MHz	6 \pm 1 Divisions			
45 MHz	6 \pm 1.8 Divisions			

Amplifier Pulse Response

Amplifier pulse response checks the aberrations, which include rise and fall times, overshoot and undershoot.

Pulse Response Checks

Equipment: Oscilloscope, PM8571A Pulse Waveform Generator

Preparation:

1. Configure the Oscilloscope follows:
Setup: As required for the test
2. Connect the 9260 input to the PM8571A output
3. Connect 9260 outputs to the oscilloscope input via 50Ω/2W Termination

Configure the PM8571A as follows:

Frequency: 1 MHz
 Waveform: Square wave
 Amplitude: 2 V
 Outputs: On

Test Procedure

1. Perform Square Wave Characteristics tests using Table 3-4.

Table 3-4, Square Wave Characteristics

Parameter Tested	Error Limits	Oscilloscope Reading	Pass	Fail
Rise/Fall Time	<10 ns			
Over/undershoot	<10%			

Auxiliary Inputs

Auxiliary inputs checks the functionality of the rear panel Inputs, which include Auxiliary Input and Auxiliary Offset.

Auxiliary Offset Checks

Equipment: Oscilloscope, PM8571A Pulse Waveform Generator, Function Generator

Preparation:

1. Configure the Oscilloscope as follows:
 - Time Base: 1 μ s/div
 - Amplitude: 5 V/div
 - Measurement: DC offset
2. Connect the 9260 input to the PM8571A output
3. Connect 9260 Channel outputs to the Oscilloscope input via 50 Ω /2W Termination
4. Configure the PM8571A as follows:
 - Amplitude: 1 V
 - Outputs: On
 - Frequency: 1MHz
5. Connect the Function Generator output to the rear panel Input
6. Configure the Function Generator as follows:
 - Wave shape: DC
 - Output: On
 - DC Amplitude: As specified in Table 3-5

Test Procedure

1. Perform Offset tests using Table 3-5.

Table 3-5, Auxiliary Offset Characteristics

Function Generator Setting	Error Limits	Oscilloscope DC offset Reading	Pass	Fail
-1V	-1V \pm 0.1V			
-0.5V	-0.5V \pm 0.05V			
0.5V	0.5V \pm 0.05V			
+1V	1V \pm 0.1V			

Auxiliary Input Checks

Equipment: Oscilloscope, PM8571A Pulse Waveform Generator

Preparation:

1. Configure the Oscilloscope as follows:
 - Time Base: 10 μ s/div
 - Amplitude: 5 V/div
 - Measurement: Amplitude Vpp
2. Connect the 9260 input CH1 to the PM8571A CH1 output and 9260 Auxiliary Input CH1 to the PM8571A CH2 output
3. Connect 9260 Channel outputs to the Oscilloscope input via 50 Ω /2W Termination
4. Configure the PM8571A as follows:
 - Outputs: Coupled
 - Frequency: 1MHz
 - Amplitude CH1: 1 V
 - Amplitude CH2: 1 V
 - CH1, CH2: On

Test Procedure

1. Measure the output amplitude, should be 10Vpp \pm 0.5V.

Adjustment Procedures

Use the following procedures to calibrate the Model 9260. The following paragraphs show how to set up the instrument for calibration and what the acceptable calibration limits are



NOTE

Initial factory adjustments require that the covers be removed from the instrument.

Amplifier Gain Adjustments

The gain adjustments ensure that amplifier gain within the specified range. Use this procedure if you suspect that the gain is an issue.

Gain Adjustments

Equipment: DMM, PM8571A AWG

Preparation:

5. Configure the DMM as follows:
Function: ACV
6. Connect the PM8571A to the DMM input
7. Configure the PM8571A as follows:
Wave Shape: Sin
Frequency: 10 kHz
Amplitude: 2Vpp
Output: On
8. Adjust the PM8571A output till DMM reading is 0.707V
9. Connect the 9260 input channel 1 to the PM8571A output
10. Connect 9260 Channel 1 output to the DMM input via 50Ω/2W Termination

Adjustment:

Adjust RV5 till output amplitude is within 7.07 V, ± 282 mV

11. Connect the 9260 input channel 2 to the PM8571A output
12. Connect 9260 Channel 1 output to the DMM input via 50Ω/2W Termination

Adjustment:

Adjust RV3 till output amplitude is within 7.07 V, ± 282 mV

Pulse Response Adjustments

The pulse response adjustments ensure that rise and fall times as well as aberrations are within the specified range. Use this procedure if you suspect that the pulse response is an issue.

Pulse Adjustments

Equipment: Oscilloscope, PM8571A AWG

Preparation:

1. Configure the PM8571A as follows:

Function:	Standard
Waveform:	Square
Frequency:	1MHz
Amplitude:	2Vpp
Output:	On

2. Connect the 9260 input channel 1 to the PM8571A output
3. Connect the 9260 output channel 1 to the oscilloscope input via 50 Ω /2W Termination.

Adjustment:

Adjust C77 for best pulse response (<10ns, 10% aberrations)

4. Connect the 9260 input channel 2 to the PM8571A output
5. Connect the 9260 output channel 2 to the oscilloscope input via 50 Ω /2W Termination.

Adjustment:

Adjust C35 for best pulse response (<10ns, 10% aberrations)

This Page was intentionally left blank

Appendix A

9260 SPECIFICATIONS

Main Input Characteristics

Characteristics	Description
No. of Channels	2 with single-ended outputs; 1 with differential outputs, 2 Summing, 2 DC offset
Connector	Front panel BNC
Impedance	50 Ω , 75 Ω or 1M Ω (factory configured)
Coupling	DC or AC
Damage Level	12Vp-p (-6V to +6V peaks)
Differential Mode Accuracy	4%

Input Auxiliary Characteristics

Characteristics	Description
Connector	Rear panel BNC
Impedance	50 Ω
Coupling	DC
Damage Level	12Vp-p (-6V to +6V peaks)

Offset Auxiliary Input Characteristics

Characteristics	Description
Connector	Rear panel BNC
Impedance	10k Ω
Coupling	DC
Damage Level	\pm 2V
Accuracy	7%

Output Characteristics

Characteristics	Description
No. of Channels	2
Coupling	DC or AC coupled
Connector	Front panel BNC
Gain (load 50 Ω)	X10, fixed (Can be ordered from factory with different gain setting. Bandwidth may change with different gain configuration)
Polarity	Normal
Impedance	0 Ω
Protection	Short Circuit to Case Ground & Thermal protection
Amplitude	34Vpp into 50 Ω
Max. Output Current	1A
Rise/Fall Time	<10ns (typ.)
Aberration	<10%

Frequency

Characteristics	Description												
Frequency Range	DC to 45MHz												
Bandwidth (-3dB) 50Ω Load	45MHz (typ.)												
Harmonics Distortion (typ.)													
	<table border="1"><thead><tr><th></th><th>10Vpp</th><th>25Vpp</th></tr></thead><tbody><tr><td>1MHz</td><td><-65dBc</td><td><-54dBc</td></tr><tr><td>10MHz</td><td><-50dBc</td><td><-45dBc</td></tr><tr><td>30MHz</td><td><-38dBc</td><td><-30dBc</td></tr></tbody></table>		10Vpp	25Vpp	1MHz	<-65dBc	<-54dBc	10MHz	<-50dBc	<-45dBc	30MHz	<-38dBc	<-30dBc
	10Vpp	25Vpp											
1MHz	<-65dBc	<-54dBc											
10MHz	<-50dBc	<-45dBc											
30MHz	<-38dBc	<-30dBc											

Mechanical, Environmental and Maintenance Specifications

Power supply

Characteristics	Description
Power Requirements	85V to 265V, 47-63Hz, <25W
Signal Ground	Grounded to case ground or Floating
EMC Certification	CE marked

Mechanical

Characteristics	Description
Physical Size	2U, half-rack size
Weight	Approximately 7Lbs

Environmental

Characteristics	Description
Temperature	
Operating	0°C to 50°C
Storage	-40°C to 70°C
Humidity	85% RH, non-condensing